





MORE THAN DEEP SKIN

BEST CHOICE IF YOUR SKIN HAVE:

- Uneven skin tone, color spots, chloasma
- Fine lines, wrinkles, looseness
- · Acne, acne marks
- Large pores and dark skin

If you have the above problem, then AURORA APL is your best choice!

GREAT FEELING YOU GET WHEN



It takes less than **3MINUTES** to treat a large back



LOVE their skin FROM HEAD TO TOE

You've seen our exceptional results. You know AURORA APL is the most versatile light device on the market with the ability to provide dermal renewal, permanent hair reduction and the treatment of pigmented lesions, vascular lesions and acne.

- >> Superior Results
- >> Speed & Precision
- >> Patient Comfort
- >> Unsurpassed Quality
- >> Cost Effectiveness
- >> User friendly&Flexible



appointments are **MAXIMIZED**



POWERFUL COOLING & Single-Degree Adjustability

>> UNSURPASSED QUALITY & ROI

» Largest Spot Size

With a large 15 x 45 mm spot size and a high repetition rate, one can treat more patients in less time with AURORA APL.



>> Uniform Beam

Precisely engineered parabolic optics deliver energy to the treatment area uniformly,producing no hot or cold spots,which results in perfectly consistent treatments.



Precision Integrated Cooling

For maximum patient comfort, AURORA APL employs an integrated, thermoelectrically cooled, sapphire crystal that keeps the treatment area consistently safe and cool with precise application dependent temperature control from $0 \, ^{\circ} \, \text{C} - 30 \, ^{\circ} \, \text{C}$.

Dual Lamp Technology

Dual lamp technology is adopted to increase the light guide area and improve the therapeutic effect greatly.




>> Cost-effective

Minimizes the operational costs and maximizes your ROI.

>> NEW Motion Technology

Delivers lower fluence at a high repetition rate for big skin area to provide patient comfort and efficacy.

>> USER FRIENDLY

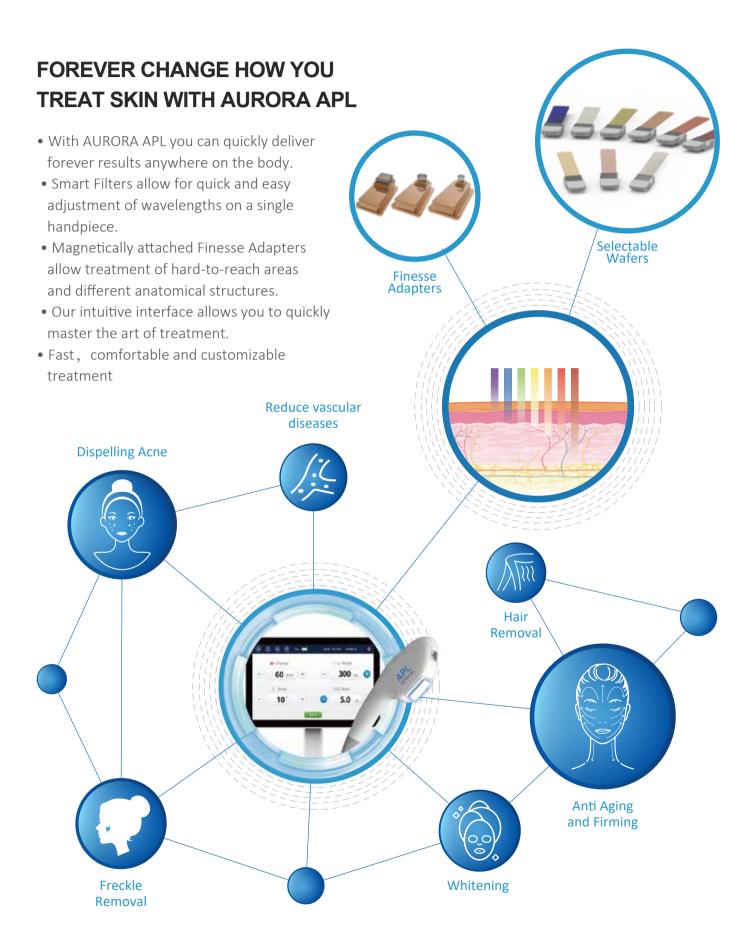

REAL PATIENTS. REAL RESULTS.

1.Ruduce Vascular **Diseases**

Hemoglobin is the main color base of blood vessels. The specific spectrum of APL photon rejuvenation can be preferentially and selectively absorbed by oxygenated hemoglobin in the blood vessels, and converted into heat energy to accumulate and heat up in the tissues. Achieve the effect of removing red blood, and

reducing vascular diseases such as acne and erythema.



Remove Red Blood Silk

2.Remove Stains

In pigmentary lesions, the selective absorption spectrum of melanin produces "implosion" or "selective pyrolysis", AURORA APL uses specific wavelengths to precisely penetrate the epidermis and dermis of the skin to destroy pigment cells and pigment particles, so as to achieve better freckle removal effect and shrink pores.



Freckle Removal. Pigment

3. Keep Skin Looking 10 Years Younger!

AURORA APL utilizes its biological stimulation to effectively improve skin conditions such as skin relaxation, wrinkles, and rough pores. increasing the activity and quantity of fibroblasts, promoting collagen synthesis, making the skin flat, smooth, and delicate.



Anti Aging and Firming

REAL PATIENTS. REAL RESULTS.

4.Forever Change How You Treat Acne

AURORA APL destroys acne-causing bacteria, reduces inflammation and redness, and makes way for clear, healthy skin. It Offers patients an alternative to expensive creams and medicine with undesired side affects. AURORA APL is fast, comfortable and non-invasive — with visible results after just one treatment, with no patient downtime.



Dispelling Acne

5.Make Skin Delicate and Whitening

AURORA APL penetrates into the deep layer of the skin through the principle of IPL, crushes the subcutaneous melanin into fine particles, and excretes it out of the body making the skin particularly delicate and whitening, and also regulating the balance of water and oil, and promoting collagen regeneration.



Whitening The Skin

6.Unwanted Hair Removal

AURORA APL uses motion technology to deliver a low fluence inmultiple pulses, protecting the skin while progressively raising the follicle temperature to the optimal level for hair removal.



Unwanted Hair Removal



YOU DESERVE **BEAUTY** YOU DESERVE **MICOWEY**

Wavelength	420 – 1400 nm
Smart Filters	420nm, 515nm, 560nm, 590nm, 640nm, 695nm, 590ST, 695ST, 800ST
Output	Up to 60 J/cm2
Pulse Width	Up to 300 ms
Repetition Rate	Up to 20 Hz
Cooling Method	Continuous thermoelectric sapphire plates
Skin Cooling	Adjustable from 0 - 30 °C in 1 °C increments
Spot Size	15 x 45 mm
Finesse Adapters	15 x 15 mm square, 11 mm & 7 mm round

MICOWEY AURORA APL





4470 OSCEOLA ST DENVER, CO 80212 USA Phone: 801.831.0172 • Email: Micowey@foxmail.com

www.micowey.com fin o







Agent Information in China

Company: Shenzhen Missdy Medical Equipment Co., Ltd

Address: 15/F, Building 10, Yuanyang Xintiandi, Shenglong Road, Longgang District, Shenzhen

Tel: 0755-32833298 Email. cadela@foxmail.com